

smiths medical
jelco®



Saf-T Wing® Blood Collection Sets

Promotes end-to-end safety for end-to-end compliance


SAF-T WING® BLOOD COLLECTION SETS

For more than 15 years, Smiths Medical has provided a wide range of innovative safety blood collection devices to help protect healthcare professionals. The Saf-T Wing® Blood Collection Set has many innovative features that help promote safety and ease-of-use.



In-vein activation

- Safety mechanism is activated while needle is in patient's vein – limiting phlebotomist's exposure to a contaminated needle
- Activated with a simple two-handed technique

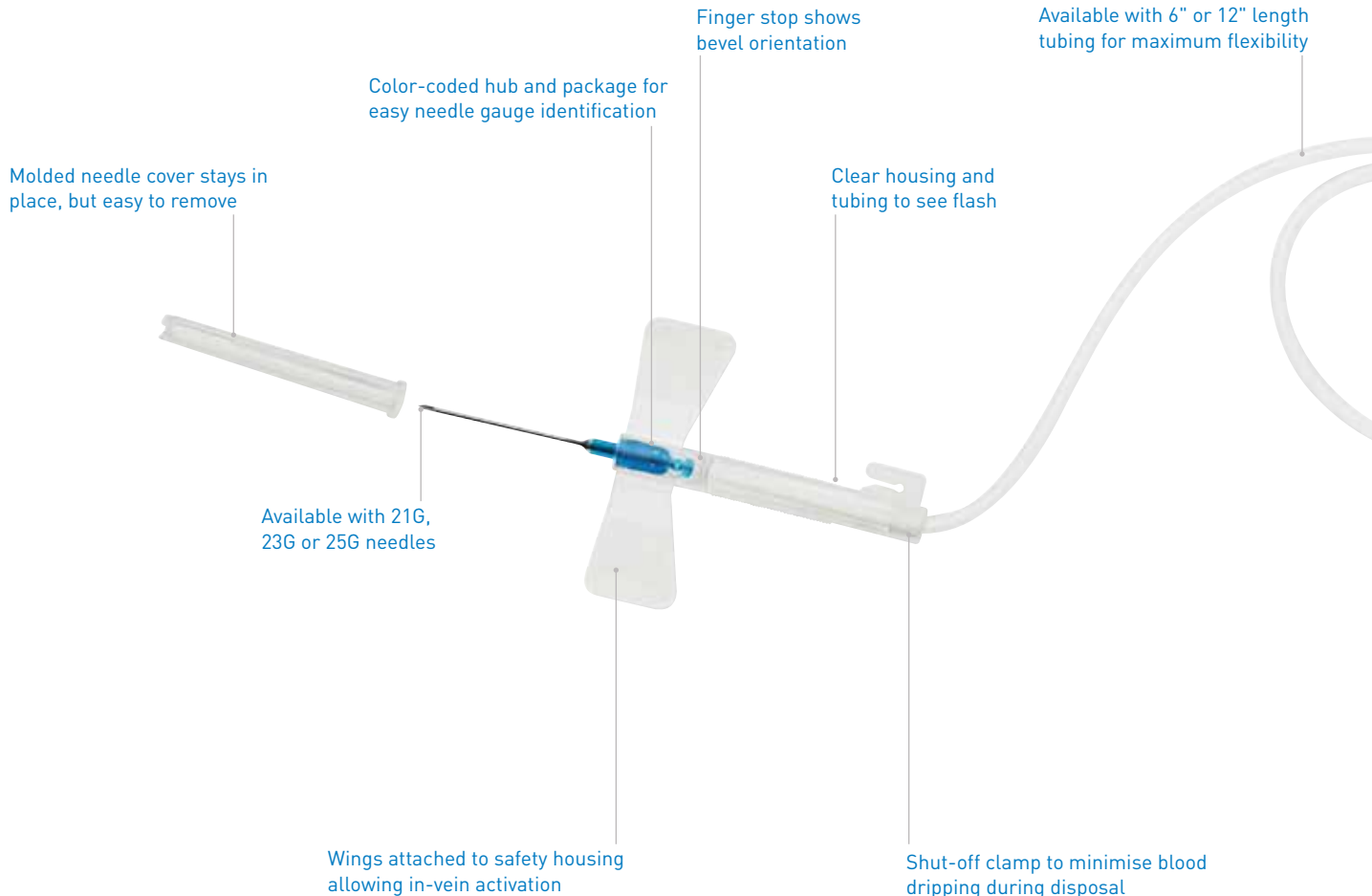



ADVANCED
FEATURES THAT
HELP PROTECT
YOU AND YOUR
PATIENTS



Slide clamp to facilitate switching
from tubing draw to syringe draw

Saf-T Pre-attached Holder® device with permanently
attached back-end needle and Luer lock to minimise
exposure and help to prevent holder re-use

Saf-T Wing® Blood Collection Sets



1

After stabilising the needle, just pull to disengage Saf-T Wing® device and activate safety feature.



2

Needle retracts and locks into place inside the safety housing when finished.



3

Wrap tubing around shut-off clamp to minimize blood dripping during disposal.

Ordering Information



Saf-T Wing® Blood Collection Sets with Tubing and Attached Holder

Product Code	Description	Color Coded	Units Per Box	Units Per Case
972106	21g x 3/4" with 6" Tubing	Green	50	200
972112	21g x 3/4" with 12" Tubing	Green	50	200
972306	23g x 3/4" with 6" Tubing	Blue	50	200
972312	23g x 3/4" with 12" Tubing	Blue	50	200
972506	25g x 3/4" with 6" Tubing	Orange	50	200
972512	25g x 3/4" with 12" Tubing	Orange	50	200



Saf-T Wing® Blood Collection Sets with Tubing and Detached Holder

982106D	21g x 3/4" Saf-T Wing® with 6" Tubing and Detached Holder	Green	50	200
982112D	21g x 3/4" Saf-T Wing® with 12" Tubing and Detached Holder	Green	50	200
982306D	23g x 3/4" Saf-T Wing® with 6" Tubing and Detached Holder	Blue	50	200
982312D	23g x 3/4" Saf-T Wing® with 12" Tubing and Detached Holder	Blue	50	200
982506D	25g x 3/4" Saf-T Wing® with 6" Tubing and Detached Holder	Orange	50	200
982512D	25g x 3/4" Saf-T Wing® with 12" Tubing and Detached Holder	Orange	50	200

All Saf-T Holder® devices incorporate a multi-sample Luer adapter



Saf-T Holder® Devices and Blood Transfer Devices

Product Code	Description	Units Per Box	Units Per Case
96000	Saf-T Holder® with Luer Adapter	50	200
96000S	Short Saf-T Holder® with Male Luer Adapter	50	200
96002	Saf-T Holder® with Female Luer Adapter	50	200
96002S	Short Saf-T Holder® with Female Luer Adapter	50	200
96004	Saf-T Holder® Blood Culture Device with Male Luer	-	50
96005	Saf-T Holder® Blood Culture Device with Female Luer	-	50
96007	Pediatric Tube Insert	10	100

All Saf-T Holder® devices incorporate a multi-sample Luer adapter

PRODUCT(S) DESCRIBED MAY NOT BE LICENSED OR AVAILABLE FOR SALE IN CANADA AND OTHER COUNTRIES

Smiths Medical International Ltd.
1500 Eureka Park, Lower Pemberton,
Ashford, Kent, TN25 4BF
Phone: +44 (0)1233 722351
Fax: +44 (0)1233 722153
www.smiths-medical.com

MKEECA-0197

smiths medical

Find your local contact information at: www.smiths-medical.com/customer-support/contact-us

Smiths Medical is part of the global technology business Smiths Group plc. Please see the Instructions for Use/Operator's Manual for a complete listing of the indications, contraindications, warnings and precautions. Saf-T Wing, Saf-T Holder and the Jetco and Smiths Medical design marks are trademarks of Smiths Medical. The symbol © indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. All other names and marks mentioned are the trademarks or service marks of their respective owners. ©2019 Smiths Medical. All rights reserved. SS195222GB-052019